

Model 2840

# Portable ceiling lift

Multiple transfer locations with one motor



SUREHANDS®





- Install rail in rooms where transfers are needed (bed, bath, toilet, shower...).
- No doorway modifications needed.
- Light, compact, only 17 ½ lbs, yet powerful: lifts 400 lbs.
- Streamlined design fits home aesthetics.
- Economical: multiple transfers with one ceiling lift.
- Use with Body Support or slings.
- 10-year design life that you can count on.
- Option: transport trolley for the Portable ceiling lift.

Specialists in innovative and reliable lift & care systems


[www.surehands.com](http://www.surehands.com)

# Portable ceiling lift

## Dimensions

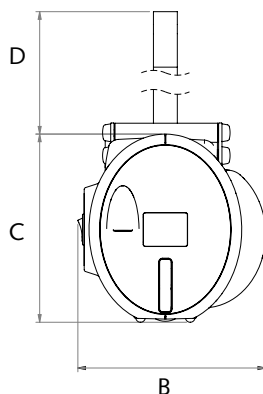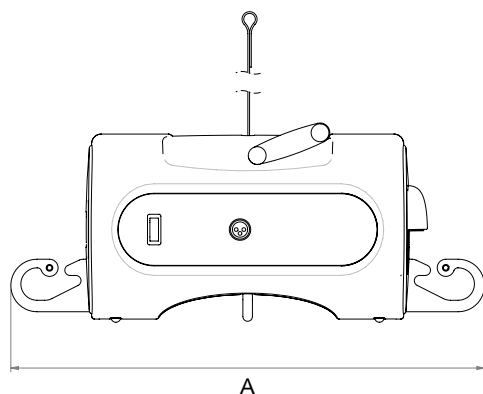

A	18.89 inch - 480 mm
B	7.8 inch - 200 mm
C	8.27 inch - 210 mm
D	61.02 inch - 1550 mm

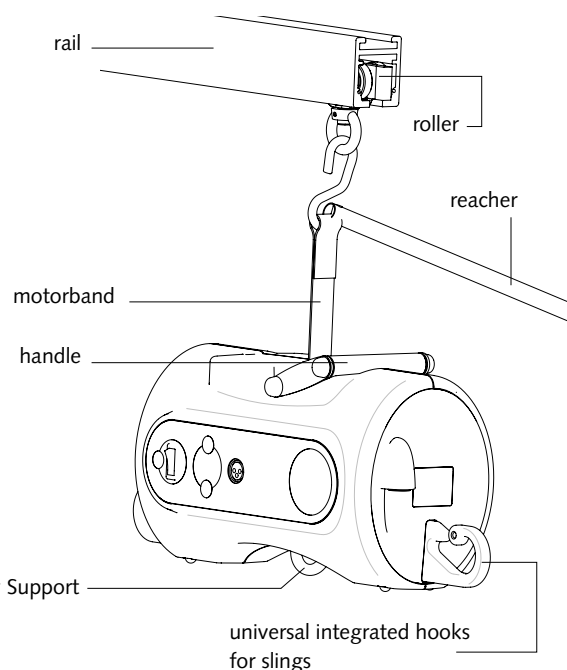


## Technical specifications

Weight .....	17.5 Lbs - 8 kg
Max. Lifting Capacity .....	400 Lbs - 182 kg
Batteries .....	2 x 12 VDC/3.2 - 3.6 Ah
Charger: .....	input 220 V AC / output 24 V DC
Number of lifts with full batteries.....	50 x with 220 Lbs - 100 kg over a distance of 19.5 inch - 50 cm
Motor cover .....	ABS
Emergency down .....	Electric
Overhead rail composition.....	Aluminum alloy, white laquered

## With Handi-slings or SureHands® Body Support

Thanks to the integrated hooks on both sides of the motor and the swivel hook placed centrally at the bottom, the Portable ceiling lift can be used with all Handi-slings or with the Body Support.

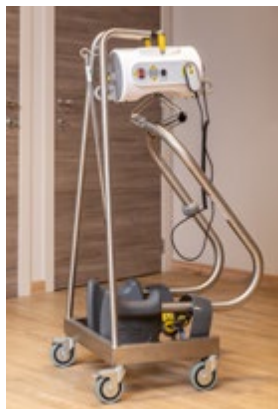

Option: Transport trolley

Local Representative

Lift & Care® is a registered trademark of  
T.F. Herceg, Inc. dba SUREHANDS LIFT & CARE SYSTEMS.

© 28-04-2025 • Handi-Move • v1.4

**SUREHANDS LIFT & CARE SYSTEMS**  
982 County Route 1, Pine Island, NY 10969  
Tel: (800) 724-5305 • Fax: (845) 258-6634 • Website: [www.surehands.com](http://www.surehands.com)