

Model 3200

# Wheelchair-to-Water™

Direct from wheelchair to pool



SUREHANDS®



- Direct transfers to and from wheelchair to water.
- Easy to move and store the lift.
- Use one lift at multiple locations: multiple pools or several areas at same pool.
- Variety of accessories make transfers easy regardless of disability.
- Designed for the swimming pool environment.


Specialists in innovative and reliable lift & care systems

[www.surehands.com](http://www.surehands.com)

## Wheelchair-to-Water™

### DIRECT TRANSFERS TO AND FROM WHEELCHAIR TO WATER.

The Wheelchair-to-Water Lift brings the user to the water in one smooth motion. This efficient lift moves the user to the pool in one transfer.

 The Wheelchair-to-Water Lift has a pivoting range of 360°.

### EASY TO MOVE AND STORE THE LIFT.

Operation is simple, with the push of a button, the spindle lowers into the floor socket at the edge of the pool.

When finished, lift the spindle out of the socket, and wheel the lift to storage.

Once removed, seal the socket opening with the plug.

### USE ONE LIFT AT MULTIPLE LOCATIONS: MULTIPLE POOLS OR SEVERAL AREAS AT SAME POOL.

The installation of additional floor sockets, and the fact that it is easily movable, allows the lift to be used in different locations

### VARIETY OF ACCESSORIES MAKE TRANSFERS EASY REGARDLESS OF DISABILITY.

Customize the transfer with the right accessory.

Based on the need choose a pool lift chair, SureHands Body Support, bathing sling or standard sling.

### DESIGNED FOR THE SWIMMING POOL ENVIRONMENT.

Made of stainless steel and powered by batteries, the mobile Wheel-chair-to Water lift is designed for safe, long-term use in the pool environment.





DISCRETE SOCKET,  
WITH PLUG TO CAP SOCKET WHEN NOT IN USE



# Wheelchair-to-Water™


## ● DIMENSIONS AND SPECIFICATIONS

| Dimensions (inch / mm) |                 | Specifications      |                       |
|------------------------|-----------------|---------------------|-----------------------|
| A                      | 72.44 / 1840.00 | Safe working load   | 300 Lbs / 135.00 kg   |
| B                      | 51.38 / 1305.00 | Batteries           | 2 x 12 Vdc; 8 Ah      |
| C                      | 21.26 / 540.00  | Charger             | separate              |
| D                      | 13.78 / 350.00  | Diameter of casters | 3.94 inch / 100.00 mm |
| E                      | 66.93 / 1700.00 | Construction        | stainless steel 316   |
| F                      | 21.65 / 550.00  | Hand control        | electrical (IP67)     |
| G                      | 19.69 / 500.00  |                     |                       |
| H                      | 12.2 / 310      |                     |                       |
| I                      | 14.57 / 370     |                     |                       |

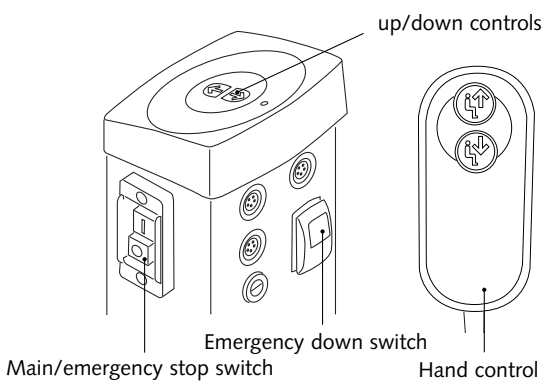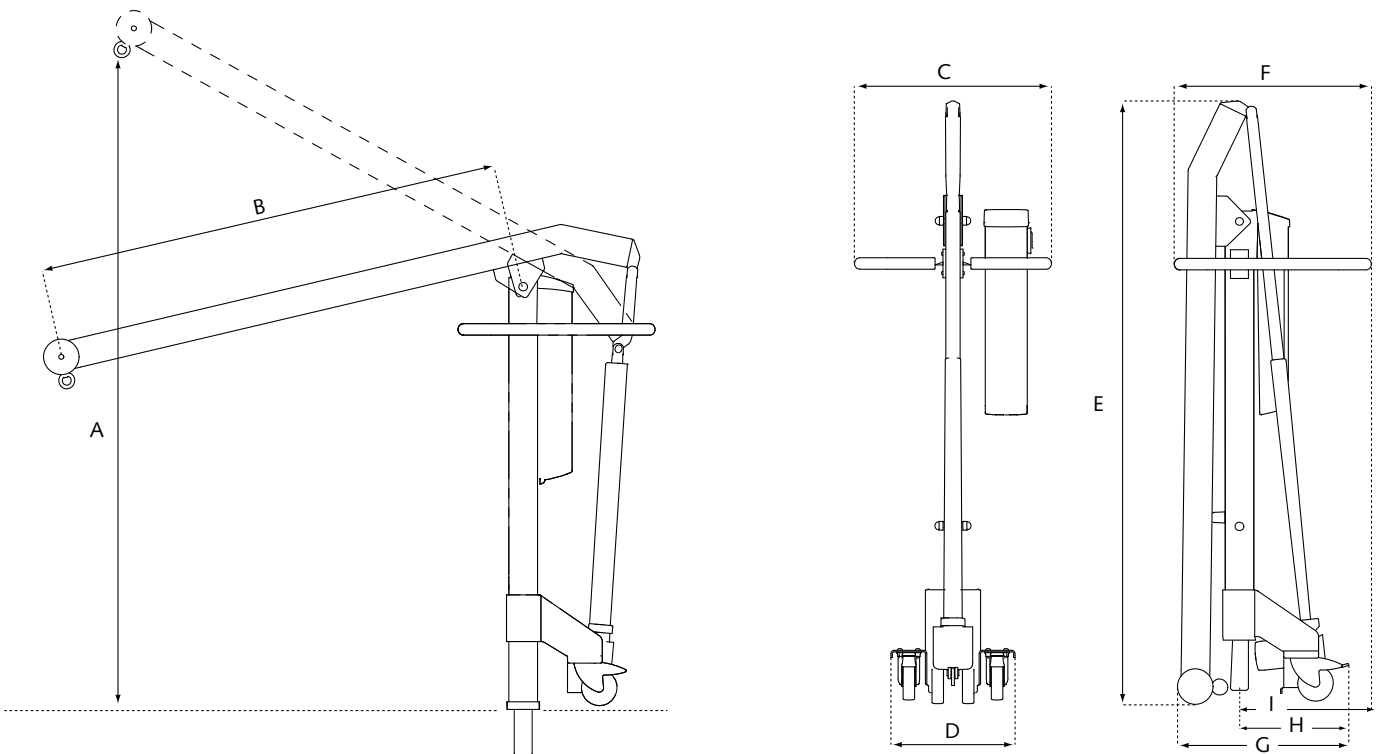



# Wheelchair-to-Water™

## ● BODY SUPPORT SYSTEMS

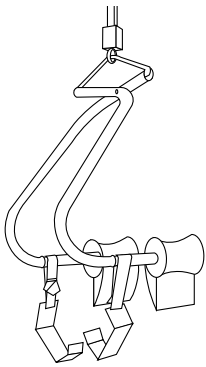



### Body Support

XS - S - M • art. 1035  
Electrolytically coated (anticorrosive)  
+ Slipproof Thigh Supports • art. 1074

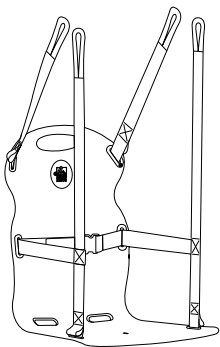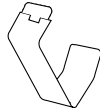

### Pool Lift Chair (PVC)



Uni-size • art. 3620  
+ Spreaderbar • art. 2110 / 2125

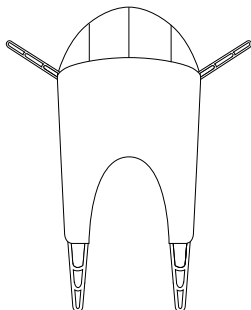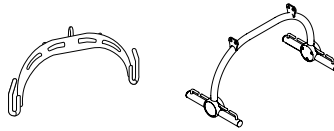

### Standard Sling


with moon shaped head support and padded thigh flaps

XS - S - M - L - XL • art. 3431  
+ Spreaderbar • art. 2110 / 2125

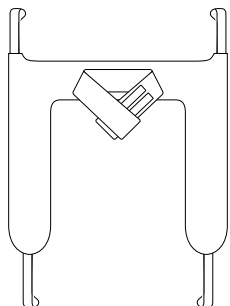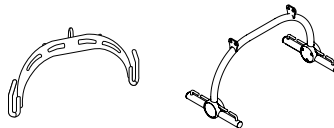

### Bathing Sling (PVC)

XS - S - M - L - XL • art. 3464  
+ Spreaderbar • art. 2110 / 2125

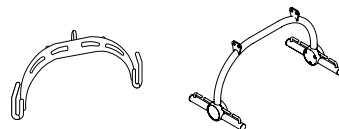

[more info](#)

Lift & Care® is a registered trademark of T.F. Herceg, Inc. dba SUREHANDS LIFT & CARE SYSTEMS.

© 31-05-2024 • Handi-Move • v1.4